

Wall-Mounted Ice and Water Dispensers

Follett's Symphony 12, 25, and 50 series of wall-mounted ice and water dispensers produce 400 pounds of ice in a 24-hour period. Between 16 and 21 inches wide, the dispensers can be mounted over a sink or anywhere space is at a premium. An integral or Satellite-fill air- or water-cooled Maestro ice-maker produces nugget ice, which rivals the cooling capacity of cube ice without the hardness that can chip teeth or the size that can cause choking. Used in ice packs, Follett's compressed nugget ice readily conforms to any shape.

All Follett wall mount dispensers are equipped with SensorSAFE infrared dispensing for enhanced sanitation; there are no levers to touch and no buttons to push. This provides less chance for cross-contamination, and allows for one-handed filling of carafes, cups, and ice bags. A companion water filter system ensures pure water quality.

► **Follett Corporation**

Circle reader service number 42
or visit www.follettice.com

Bar Code Label Printer



Seiko Instruments' SLP 450 thermal label printer can sit on any desktop surface and create high-quality labels, incorporating text, hi-res graphics, and bar codes in a wide variety of sizes. The SLP 450 can also produce 2D bar codes.

► **Seiko Instruments**

Circle reader service number 29
or visit www.seikoinstruments.com

Syringe to Protect Workers from Needle-Stick Injuries

Inviro Medical launches the InviroSNAP! Safety Syringe. With hand and fingers safely positioned behind the needle, the clinician administers the dose, pulls the plunger back — causing the needle to retract inside



the barrel — and snaps off the plunger, permanently disabling the device. The InviroSNAP! Safety Syringe delivers an accurate dose with minimal medication waste. In addition, because the needle retracts into the barrel, the syringe can decrease the risk of injury after the injection, improve infection control, and lower related disposal costs.

► **Inviro Medical**

Circle reader service number 30
or visit www.inviromedical.com

Long-Term Care EMAR Solution

Artromick launches a new EMAR (electronic medication administration record) solution specifically designed for long-term care. The Web-based e-charting solution offers a user-friendly interface and the ability to access patient information on a 24-7 basis.



► **Artromick International**

Circle reader service number 31
or visit www.artromick.com

Upgraded Compounding Software

New software for Grifols' Gri-fill System improves the device's automated compounding process for sterile IV and epidural preparations. The system is used for batch preparations such as PCA, epidurals, electrolytes, and anti-infectives, as well as custom preparations like oncology drugs. The new software allows you to automatically perform and document accuracy checks, and has new features to minimize drug waste and an optional feature for bar coding IV labels.



► **Grifols**

Circle reader service number 32
or visit www.grifolsusa.com

Narcotic Control Cabinets

EPS's line of narcotic control cabinets are designed to hold appropriate quantities of controlled substances in limited spaces. They can be placed wherever such products need to be accessible but secure, and can be mounted on walls. Some models are equipped with double doors, and all doors use separate keys for additional security. The keys cannot be removed when the lock is in the open position. Manufactured from 22-gauge stainless steel with aluminum shelves, EPS narcotic control cabinets are sturdy, durable, and ideal for narcotic security when space is at a premium.



► **EPS, Inc.**

Circle reader service number 34
or visit www.medi-dose.com

Pharmaceutical Waste Containers

Kendall's PharmaSafety Pharmaceutical Waste Containers for the disposal of non-hazardous pharmaceutical waste are blue and white in color, making them distinguishable from sharps, chemo, and other disposal containers. Containers range from 2 to 18 gallons in size and include tamper-evident labels. Leak-resistant gasketing and absorbent pads help contain liquid contents. In-room wall enclosures, wire brackets, table top holders, and locking floor carts are available.

Kendall's RCRA Hazardous Waste Containers can be used to segregate hazardous waste in accordance with EPA requirements. Containers are black in color, range from 2 to 18 gallons in size, and meet all EPA labeling requirements. Leak-resistant



gasketing and absorbent pads help contain liquid contents. Locking floor carts are available. Kendall's waste containers are available on the following GPO contracts: Amerinet, Consorta, Healthtrust, Med-Assets, MSA, Novation, and Premier.

- Kendall
Circle reader service number 33
or visit www.kendallhq.com

2D Bar Coded, Unit Dose Packaging

Medical Packaging Inc.'s Bar-Code Plus Solutions are now capable of printing 2-D RSS Symbology, which can store 100 times more information and is 30 times smaller than a linear, one-dimensional barcode. Choose from a variety of 2-D barcode formats, including Data Matrix, PDF417, and RSS Limited. All formats can be customized and typically include NDC number, lot number, and expiration date.




formats can be customized and typically include NDC number, lot number, and expiration date.

- Medical Packaging Inc.
Circle reader service number 35
or visit www.medpak.com

Tamper-Evident Safety Seals

EPS now offers a bright red Bio-Hazard Seal for securing hazardous waste. New Pull-Thru Seals for crash carts and emergency drug boxes

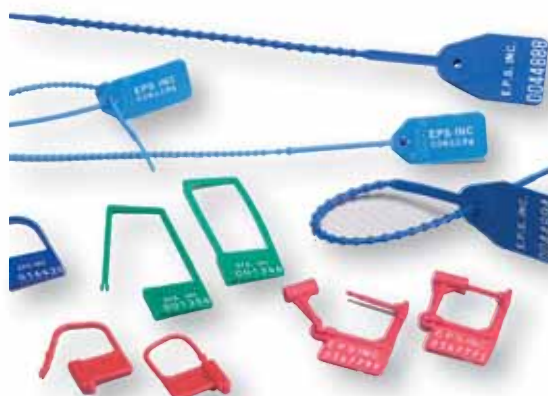

are also available, and a variety of Plastic Padlock Seals include both long and short hasps. All Plastic Padlock Seals are available either numbered or unnumbered.

- EPS, Inc.
Circle reader service number 36
or visit www.medi-dose.com

Bar Coded Wristbands

Precision Dynamics introduces Soft CompuBand patient wristbands, made of a soft material for increased patient comfort and featuring a tamper-proof, single-post snap for added ease-of-use. Soft CompuBands support text, linear and two-dimensional bar codes, patient photos, and other graphics. The wristbands are available in adult and pediatric sizes in up to six colors.



sional bar codes, patient photos, and other graphics. The wristbands are available in adult and pediatric sizes in up to six colors.

- Precision Dynamics Corporation
Circle reader service number 38
or visit www.pdc-healthcare.com

Bar Code Scanner

Hand Held Products introduces the 3800g linear bar code imager. The new scanner offers improved depth of field, faster bar code reading, and a higher tolerance to motion. The device's ergonomic design can reduce user fatigue and enhance productivity. With no internal moving parts, the 3800g comes with popular connectivity interfaces on-board and is backed by a five-year warranty.



- Hand Held Products
Circle reader service number 40
or visit www.handheld.com

Compounding Aseptic Isolator

Germfree's VersaFlow Compounding Aseptic Isolator (CAI) is specifically designed and engineered for low- to medium-risk compounding. A deep work area accommodates automated compounding devices. The stainless-steel VersaFlow can be converted to operate as a CAI or as a laminar airflow workstation (LAFW), for use inside a cleanroom.

- Germfree Laboratories
Circle reader service number 41
or visit www.germfree.com

Compounding Documentation System

The VALITEQ Documentation Management System (DMS) version 2.0 provides pharmacy supervisors with detailed records and reports on the personnel, environmental, and facility status of their compounding programs, leading to the level of process control mandated by USP <797>. The Java-based DMS v2.0 provides audits of your compounding processes, allowing you to detect and correct process parameters that may be trending out-of-limits, such as personnel competency, performance, and monitoring; environmental contamination monitoring and sanitizing; and facility engineering control (HVAC and LAFW) operation, service, and maintenance.

- Lab Safety Corporation
Circle reader service number 42
or visit www.valiteq.com

Anti-Microbial Additive for Pneumatic Tube System Carriers

Swisslog introduces BacStop, an optional anti-microbial additive for its standard TransLogic and Eco-Seal leak-resistant pneumatic tube system carriers. BacStop is a silver ion anti-microbial agent molded within the body of the pneumatic tube carrier. BacStop has been confirmed effective against more than 35 different bacteria, yeast, and fungi.

- Swisslog
Circle reader service number 64
or visit www.swisslog.com

Pharmacy Refrigerators

Helmer Pharmacy Refrigerators are available in both single- and double-door models and offer storage capacities from 20 to 56 cubic feet. The units feature a bacteria-resistant-powder-coated interior and exterior, dual-pane glass doors with key locks,

and solid stainless-steel drawers. Helmer Pharmacy Refrigerators are constructed to pass through standard height doorways and can be configured with a combination of drawers; zinc-plated, powder-coated wire shelves; and roll-out wire baskets. Helmer's i.Series line features the i.Center Integrated Monitoring System, which constantly monitors the chamber temperature and other critical refrigerator functions to safeguard the products inside and provide information to users.

- Helmer Inc.
Circle reader service number 63
or visit www.helmerinc.com

Unit Dose Packaging Services

Baxa and Unit Dose Solutions (UDS) have signed an exclusive co-marketing agreement for the promotion of UDS' services for hospital pharmacy drug repackaging. Baxa will market and sell outsourced repackaging services on behalf of UDS, and provide oral and topical dispensers and accessories for its operations. Using current Good Manufacturing Practices, Unit Dose Solutions can repackage medications in a unit dose, bar-coded format. UDS works with the customer's existing supply chain, using their existing drug formulary, and coordinating with their drug wholesaler and their internal storage and distribution system. UDS can also provide individual product labeling for a variety of bedside verification systems and bar code scanners.

- Baxa Corporation
Circle reader service number 62
or visit www.baxa.com

Compounding Aseptic Isolator

Baker's new SS800 compounding aseptic isolator has two double-glove work chambers connected by an internal pass through, allowing product hand-off between two technicians in multi-step compounding activities. In addition, the internal pass-through doors may be closed, allowing each chamber to operate as a stand-alone isolator. Like Baker's other SterilSHIELD models, the SS800 utilizes laminar airflow within each of the main work chambers and pass-through chambers for rapid removal of contaminants, and features a sloped front panel for user comfort. The outer pass-through door can allow for the efficient transfer of materials.

- The Baker Company
Circle reader service number 61
or visit www.bakerco.com



Handwash Kits

Glitterbug Handwash Show kits from Brevis Corporation can aid in your facility's hand hygiene education efforts. Using the fluorescent GlitterBug Potion and a UVA lamp to illuminate

present microorganisms, your staff can see how well their hands are washed. GlitterBug Powder is used to demonstrate how contamination and germs can travel. The starter kit includes a GlowBar lamp with four D-Cell batteries, an AC adapter for plug-in convenience, a hand wash

DVD, GlitterBug Potion (240ml), GlitterBug Powder (68gm), an instruction manual, a "How to Wash" poster, and two hand brushes.

► **Brevis Corporation**

Circle reader service number 125
or visit www.brevis.com

Pharmaceutical Waste Stream Analysis

PharmEcology's PharmE Formulary Analysis allows you to gain a clear view of your pharmaceutical waste stream with an immediate and complete analysis of your entire inpatient formulary. PharmEcology can provide you with complete pharmaceutical waste stream management recommendations, identifying all federally hazardous and PharmE Hazardous waste

present in your formulary and inventory. PharmE Hazardous waste classifies products which, while not technically hazardous, should be treated as hazardous waste due to their potential to cause harm. Using PharmE Formulary Analysis, you can demonstrate active compliance efforts by identifying all hazardous pharmaceutical waste within 10 working days.

► **PharmEcology Associates**
Circle reader service number 124
or visit www.pharmecology.com



Investigational Drug Service Management System

WebIDS by the McCreadie Group automates many of the manual processes currently performed in an investigational drug service (IDS).



WebIDS manages the research protocols and handles the dispensing, inventory, and record-keeping needs of an IDS. Fully automated billing, dispensing labels, and inventory-accountability logs reduce burdensome manual tasks. These features, along with IRB tracking, automated logging of records changes, contact management, and reporting, help streamline the daily operations of an IDS.

► **McCreadie Group, Inc.**

Circle reader service number 123
or visit www.mccreadiegroup.com

Robotic IV Automation

The RIVA system automates IV preparation in compliance with USP <797> and other regula-

tions. RIVA has been designed to safely prepare non-hazardous drugs (including pediatric doses), as well as chemotherapy and other hazardous drugs. Integrated with the pharmacy

information system, RIVA fulfills pharmacist-approved medication orders in a fully enclosed ISO Class 5 environment with continuous air particle monitoring. RIVA prepares syringes from 1 to 60 mL and bags from 50 to 1000 mL, and reconstitutes drug vials from 1 mL to 100 mL. A permanent electronic record is created with every preparation, and customized labels are affixed to preparations inside the device. RIVA can prepare both batch and patient-specific orders, and production can be stopped at any time for STAT orders. The robot's estimated throughput is 800 doses in an eight-hour shift.

► **Intelligent Hospitals Ltd.**

Circle reader service number 122
or visit www.rivasystem.com



New Intelligent Infusion Pump

Hospira unveils the Symbiq infusion system, an "intelligent" pump that incorporates clinician-friendly features to decrease errors and improve



workflow. Hospira's MedNet software is integrated with the pump, offering IV administration safeguards and sophisticated reporting capabilities. The pump's LCD display features simple layouts, large fonts, color prompts, and touch-activated buttons to help minimize errors, and is visible from distances of up to 12 feet. A quick-release pole clamp enables clinicians to rapidly remove or attach the system to an IV pole or bed rail. An automatic medication-loading slot allows for additional ease of use and enhanced free-flow protection.

► **Hospira, Inc.**

Circle reader service number 120
or visit www.hospira.com

Luer Lock Closures

EPS announces the addition of another color to its growing line of Combi-Cap Male/Female Luer Lock closures. Combi-Cap closures offer security and versatility when capping syringes, catheters, and other hypodermic devices. A new blue cap has been added to complement the existing red, white, and orange caps. Combi-Cap closures allow pre-filled syringes to stand upright.

One side can act as a syringe cap, and the other as a catheter cut-off. Individual packaging can prevent compromises in sterile handling.

► **EPS, Inc.**

Circle reader service number 121
or visit www.medidose.com

